


### Mechanical Lung (Rear View)

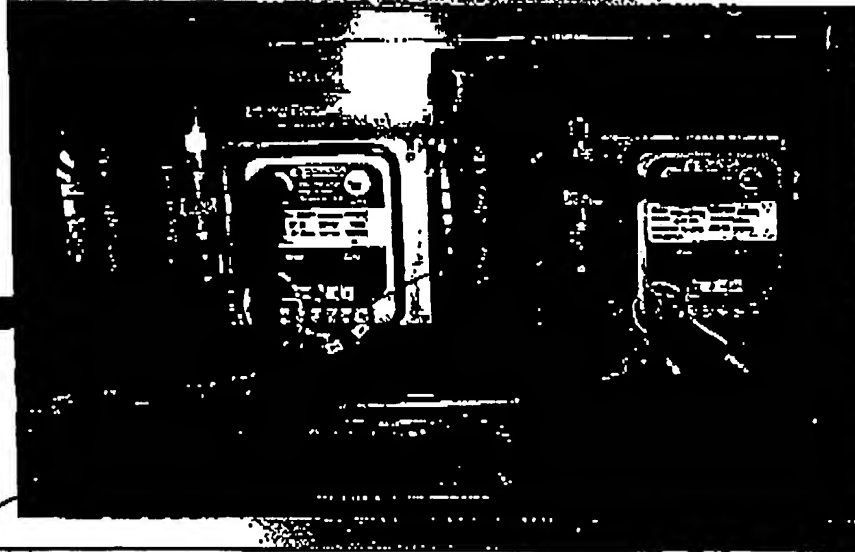


### Inhalation

Air Enters the  
Nose and Mouth

The Chest  
Elastically Expands

Diaphragm Is  
Forcefully Lowered

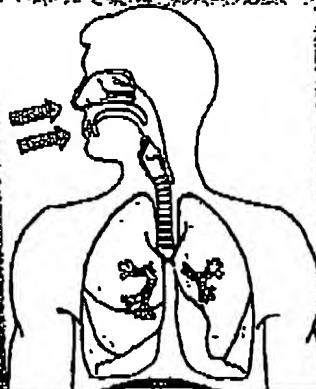

Resistance Is  
Offered in the  
Trachea and


FIG. 1

## Exhalation

Air Exits the  
Nose and Mouth

The Chest Tends  
to Be Un-Deformed  
Position

The Diaphragm  
is Relaxed (Passive)

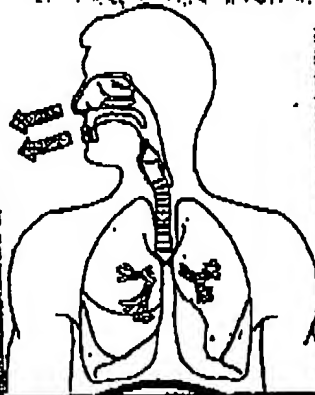


Resistance is  
Created in the  
Lungs and

FIG. 2

## Respiration Control

Diaphragm

Linear  
Actuator

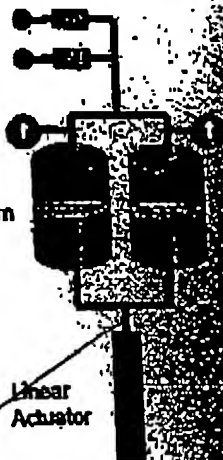

UNITED STATES GOVERNMENT  
 OFFICE OF THE SECRETARY OF DEFENSE  
 WASHINGTON, D. C. 20301-4000

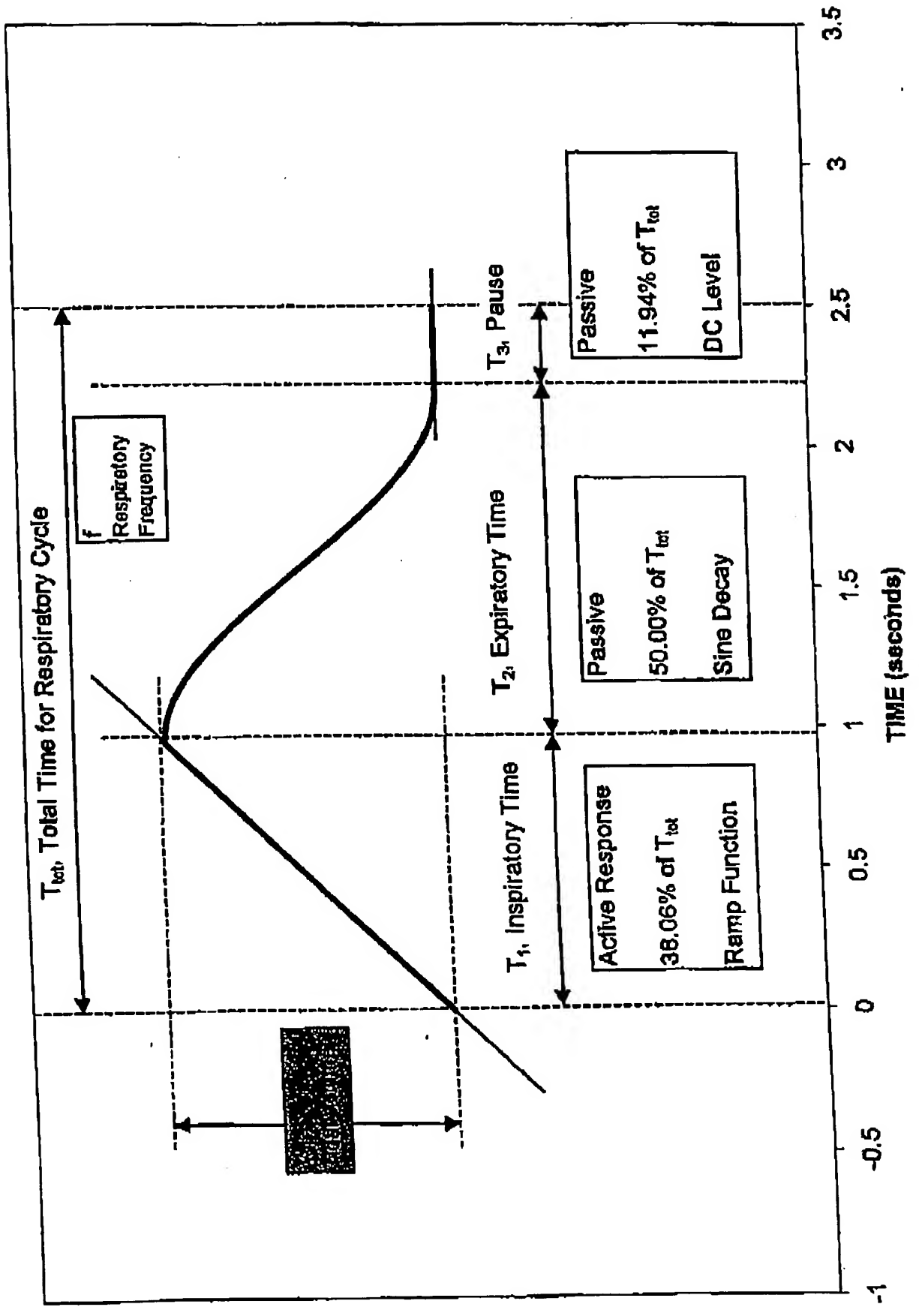


FIG. 3

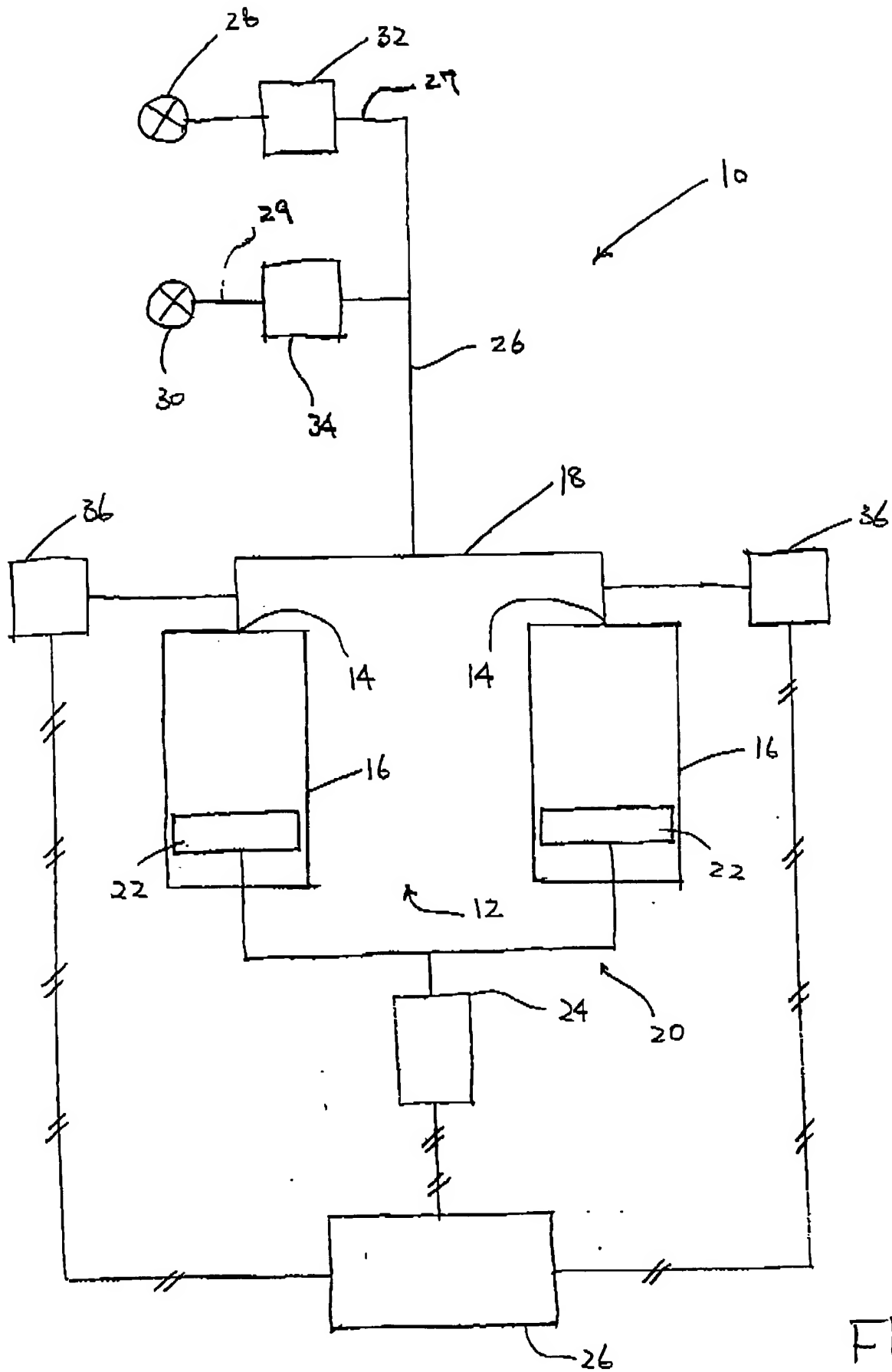
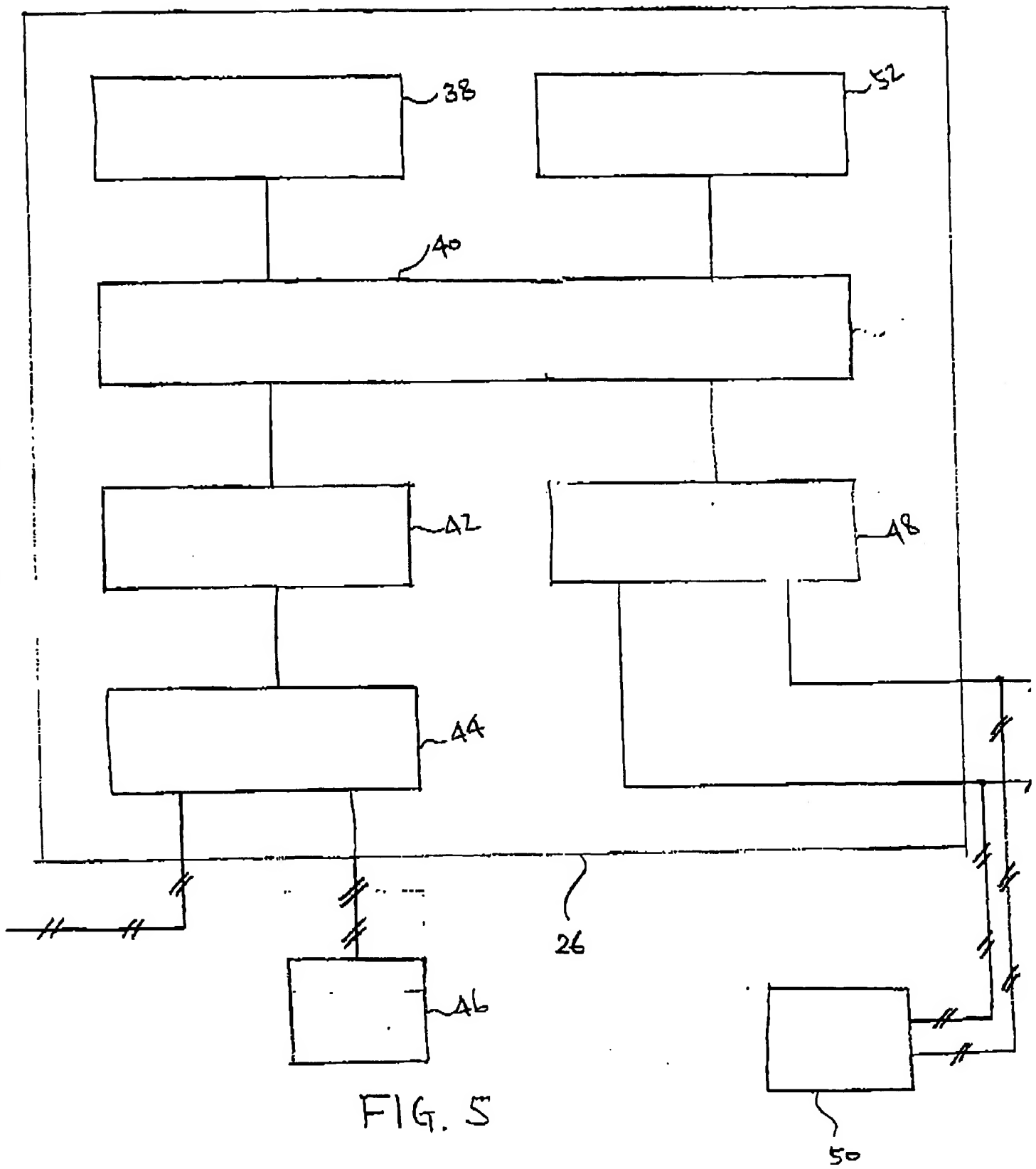


FIG. 4



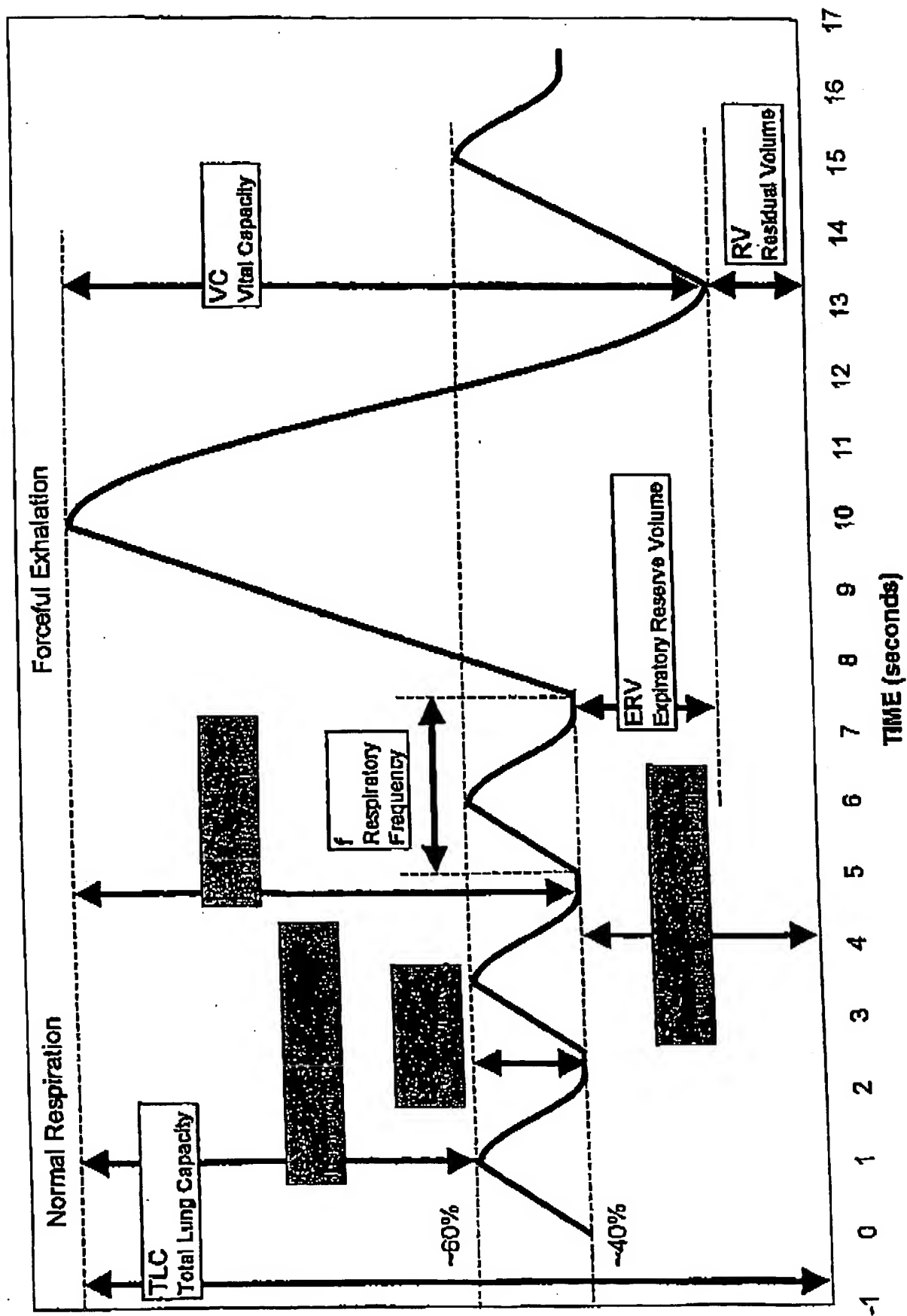


FIG. 6

NOTE: The values in this table are for a normal adult male. For other subjects, the values should be adjusted accordingly.

AGE (years)												
	1	2	3	4	5	6	7	8	9	10	11	12
Lung Capacity												
TLC (ml)	833	823	1215	1340	1467	1802	2138	2473	2788	3123	3448	3773
Total Lung Capacity												
VC (ml)	475	683	910	1005	1100	1362	1603	1855	2089	2343	2586	2830
Vital Capacity												
RV (ml)	158	231	303	335	367	451	534	618	700	781	862	843
Residual Volume												
TLC = VC + RV												
RV = 0.25 * TLC												
Normal Respiration												
Functional Residual Capacity												
FRC (ml)	263	388	532	668	860	831	1003	1174	1344	1515	1685	1855
VT (ml)	78	95	112	121	130	147	163	180	200	220	240	260
Tidal Volume												
IC (ml)	370	528	681	744	807	871	1135	1298	1454	1608	1764	1918
Inspiratory Capacity												
IRV (ml)	282	431	589	623	677	824	872	1119	1254	1388	1524	1658
Inspiratory Reserve Volume												
Frequency (cycles/min)	24	23	22	21	20	19	18	18	18	17	17	16
TLC = FRC + VT												
TLC = FRC + VT + IRV												
IC = VT + RV												
FRC = 0.30 * TLC (upright)												
FRC = 0.40 * TLC (supine)												
Foreaful Exhalation												
ERV (ml)	105	167	228	281	283	361	468	558	645	734	823	912
Expiratory Reserve Volume												
Frequency (cycles/min)												
TLC = IC + ERV + RV												
VC = ERV + IC												

Note: Volume = the volume of 1 lung, Tot Vol = the volume of both lungs

Enter values at "red" tent only.

Position	Volume (ml)	Tot Volume (ml)
Baseline (0.00-inches, 0.0 volts)	0.00	0.00
Fully Extended (6.00-inches, 5.0 volts)	8.00	8.00
V <sub>max</sub>		183

mAV													
Normal Respiration Cycle													
Total Time (1 Respiration Cycle)	T <sub>tot</sub> (sec)	2.508	2.608	2.717	2.857	3.000	3.103	3.214	3.333	3.428	3.528	3.636	3.744
Inspiratory Time (T <sub>i</sub> ) Ramp, 38.05%	T <sub>i</sub> (sec)	0.952	0.993	1.038	1.087	1.142	1.181	1.223	1.268	1.305	1.343	1.384	1.423
Expiratory Time (T <sub>e</sub> -T <sub>i</sub> ) Sine Decay, 50%	T <sub>e</sub> (sec)	1.250	1.304	1.384	1.429	1.500	1.552	1.607	1.687	1.714	1.785	1.818	1.943
Pause (T <sub>e</sub> -T <sub>i</sub> ) DC Level, 11.94%	T <sub>p</sub> (sec)	0.299	0.311	0.328	0.341	0.358	0.371	0.384	0.398	0.408	0.421	0.434	0.459
Programming Settings:													
RAMP	RS (V)	0.2814	0.6489	1.0164	1.1913	1.3881	1.8342	2.3024	2.7705	3.2357	3.7008	4.1680	4.6352
Start	RE (V)	0.4845	0.8085	1.3224	1.5218	1.7213	2.2350	2.7488	3.2623	3.7621	4.3018	4.8217	5.3416
End													
SINE	SA (V)	0.1088	0.1288	0.1530	0.1653	0.1778	0.2004	0.2231	0.2458	0.2732	0.3005	0.3278	0.3552
Amplitude	SD (V)	0.3880	0.7787	1.1884	1.3566	1.5437	2.0348	2.5255	3.0184	3.5089	4.0014	4.4939	4.9864
Offset	SF (Hz)	0.4000	0.3833	0.3667	0.3500	0.3333	0.3222	0.3111	0.3000	0.2889	0.2833	0.2750	0.2687
Freq (Hz)	SP (deg)	80	90	90	90	90	90	90	90	90	90	90	90
phase													
DC LEVEL	DO (V)	0.2814	0.6489	1.0164	1.1913	1.3881	1.8342	2.3024	2.7705	3.2357	3.7008	4.1680	4.6352
Offset													

FIG. 7

## Manikin Connections

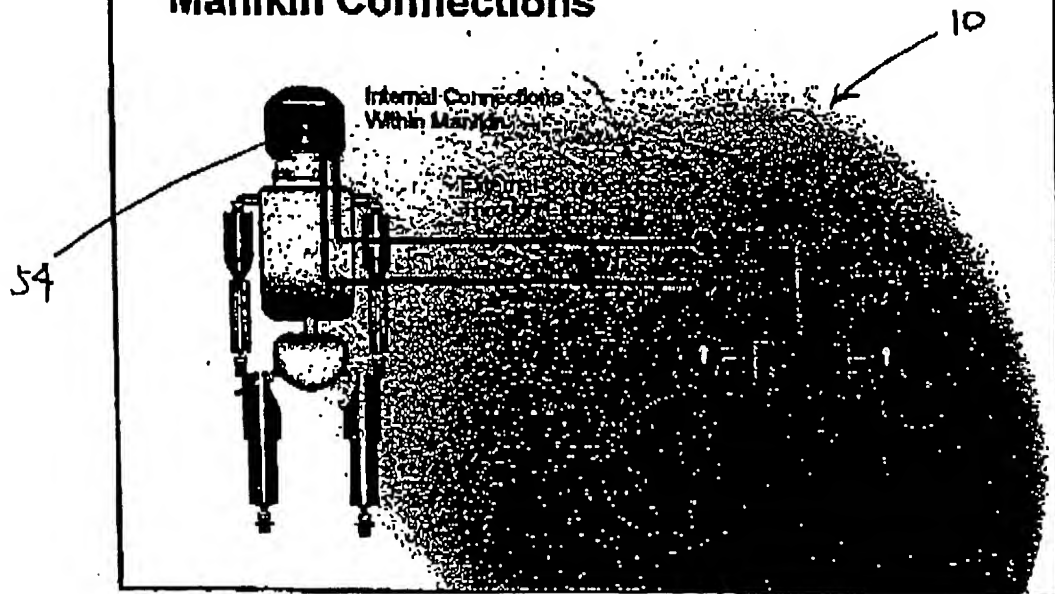


FIG. 8

## Manikin Connections

